




563 Ionization Gauge Tube

INSTRUCTION MANUAL

Manual No. 699905563
Revision B
July 2002

563 Ionization Gauge Tube

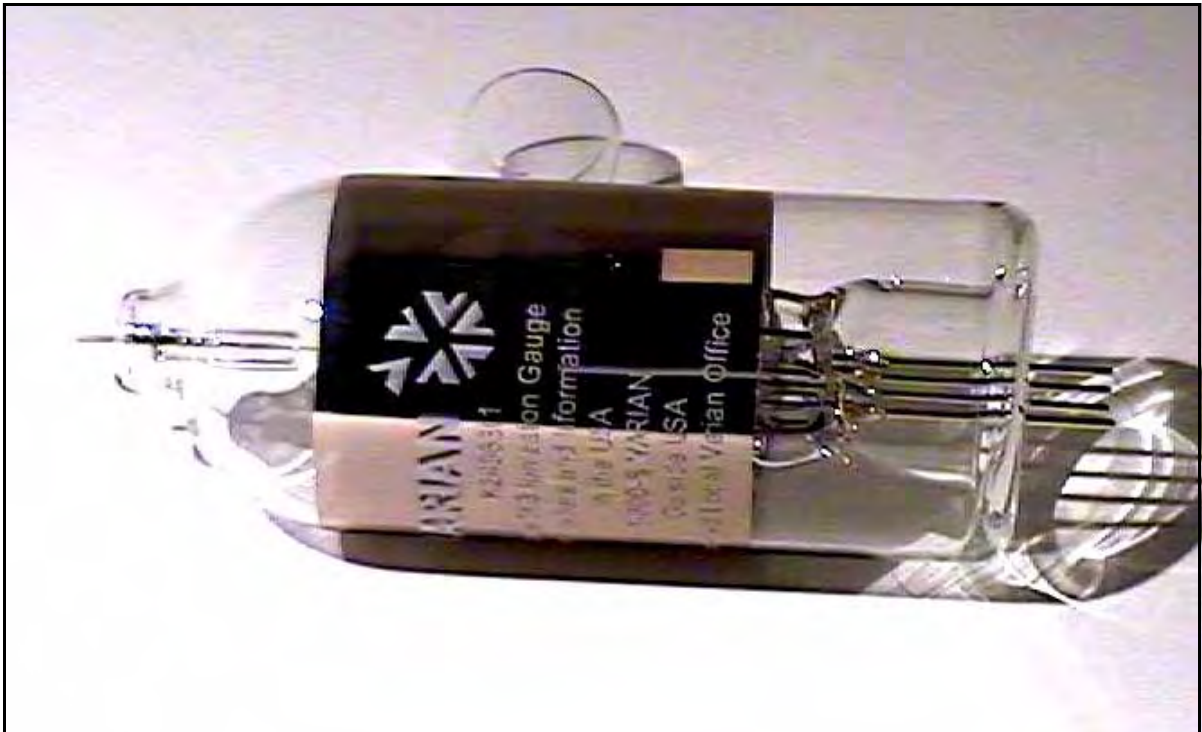

Warranty

Products manufactured by Seller are warranted against defects in materials and workmanship for twelve (12) months from date of shipment thereof to Customer, and Seller's liability under valid warranty claims is limited, at the option of Seller, to repair, to replace, or refund of an equitable portion of the purchase price of the Product. Items expendable in normal use are not covered by this warranty. All warranty replacement or repair of parts shall be limited to equipment malfunctions which, in the sole opinion of Seller, are due or traceable to defects in original materials or workmanship. All obligations of Seller under this warranty shall cease in the event of abuse, accident, alteration, misuse, or neglect of the equipment. In-warranty repaired or replaced parts are warranted only for the remaining unexpired portion of the original warranty period applicable to the repaired or replaced parts. After expiration of the applicable warranty period, Customer shall be charged at the then current prices for parts, labor, and transportation.

Reasonable care must be used to avoid hazards. Seller expressly disclaims responsibility for loss or damage caused by use of its Products other than in accordance with proper operating procedures.

Except as stated herein, Seller makes no warranty, express or implied (either in fact or by operation of law), statutory or otherwise; and, except as stated herein, Seller shall have no liability under any warranty, express or implied (either in fact or by operation of law), statutory or otherwise. Statements made by any person, including representatives of Seller, which are inconsistent or in conflict with the terms of this warranty shall not be binding upon Seller unless reduced to writing and approved by an officer of Seller.

Warranty Replacement and Adjustment

All claims under warranty must be made promptly after occurrence of circumstances giving rise thereto, and must be received within the applicable warranty period by Seller or its authorized representative. Such claims should include the Product serial number, the date of shipment, and a full description of the circumstances giving rise to the claim. Before any Products are returned for repair and/or adjustment, written authorization from Seller or its authorized representative for the return and instructions as to how and where these Products should be returned must be obtained. Any Product returned to Seller for examination shall be prepaid via the means of transportation indicated as acceptable by Seller. Seller reserves the right to reject any warranty claim not promptly reported and any warranty claim on any item that has been altered or has been returned by non-acceptable means of transportation. When any Product is returned for examination and inspection, or for any other reason, Customer shall be responsible for all damage resulting from improper packing or handling, and for loss in transit, notwithstanding any defect or non-conformity in the Product. In all cases, Seller has the sole responsibility for determining the cause and nature of failure, and Seller's determination with regard thereto shall be final.

If it is found that Seller's Product has been returned without cause and is still serviceable, Customer will be notified and the Product returned at its expense; in addition, a charge for testing and examination may be made on Products so returned.

3/1/00

Hazard and Safety Information

This manual uses the following standard safety protocols:

WARNING



The warning messages are for attracting the attention of the operator to a particular procedure or practice which, if not followed correctly, could lead to serious injury.

CAUTION



The caution messages are displayed before procedures, which if not followed, could cause damage to the equipment.

NOTE



The notes contain important information.

Operators and service personnel must be aware of all hazards associated with this equipment. They must know how to recognize hazardous and potentially hazardous conditions, and know how to avoid them. The consequences of unskilled, improper, or careless operation of the equipment can be serious. This product must only be operated and maintained by trained personnel. Every operator or service person must read and thoroughly understand operation/maintenance manuals and any additional information provided by Varian Vacuum Technologies. All warning and cautions should be read carefully and strictly observed. Consult local, state, and national agencies regarding specific requirements and regulations. Address any safety, operation, and/or maintenance questions to your nearest Varian Vacuum Technologies office.

Use with Combustibles and Mixtures

WARNING



As with all ionization gauges, this device is not intrinsically safe. Exercise extreme care when using this vacuum gauge while pumping or backfilling a system or in any other system condition which contains combustible gases or mixtures. The filament of a hot filament ion gauge and the high voltage discharge of a cold cathode gauge can be ignition sources.

When such a gas or mixture is present, do not turn on any such vacuum gauge.

Failure to follow this instruction could result in serious injury to personnel and damage to equipment.

Vacuum Equipment Cleanliness

Cleanliness is vital when servicing any vacuum equipment.

CAUTION



Do not use silicone oil or silicone grease.

Use powder-free butyl or polycarbonate gloves to prevent skin oils from getting on vacuum surfaces.

Do not clean any aluminum parts with Alconox[®]. Alconox is not compatible with aluminum and will cause damage.

NOTE



Normally, it is unnecessary to use vacuum grease. However, if it must be used, do not use silicone types, and use it sparingly. Apiezon[®] L grease is recommended (Varian Part Number 695400004).

Contacting Varian Vacuum Technologies

In the United States, you can contact Varian Vacuum Technologies Customer Service at 1-800-8VARIAN.

Internet users:

- ❑ Send email to Customer Service & Technical Support at vpl.customer.support@varianinc.com
- ❑ Visit our web site at www.varianinc.com/vacuum
- ❑ Order on line at www.evarian.com

See the back cover of this manual for a listing of our sales and service offices.

Introduction

The Varian 563¹ Ionization Gauge Tube is a broad-range Bayard-Alpert device for measuring vacuum in a wide range of vacuum pressures. Careful adherence to manufacturing techniques and standards insures close agreement between gauges. Initial calibration against closely checked McLeod gauge standards insures that this calibration is as well founded as present state of the art will permit.²

Unique with its high-pressure capability and low X-ray limit, the 563 tube is usable to 1×10^{-1} Torr in argon and 6×10^{-2} Torr in air. This is accomplished by the smaller geometry of the tube elements and the use of a conductive platinum coating on the inside of the glass envelope which improves gauge tube performance.

In high and ultra-high vacuum ranges, a platinum shield prevents external electrical fields from affecting the low ion current. In higher pressure ranges where the number of ions produced is much greater, the platinum conductive coating drains off the static charge which can build up on the inside of the glass tube. Thus, the 563 tube gives superior performance throughout the entire range of the tube.

WARNING



Check that your Ion Gauge Controller and vacuum system are separately grounded to a common ground.

Placement of a ground wire between the vacuum chamber and the controller chassis is not safe; large continuous currents could flow through it.

Personnel could be killed by high voltages (160 to 1000 volts) which may be present in an improperly grounded system).

Make absolutely sure that your vacuum system is grounded as shown in the following schematic diagram; test the system ground to be sure that it is complete and capable of supporting at least 10 Amps.

1. US Patent 3153744.

2. See: "A Detailed Examination of the Principles of Ion Calibration", W.B. Nottingham and F.L. Tomey, Jr. Vacuum Symposium Transactions, 1960, p.117.

563 Ionization Gauge Tube

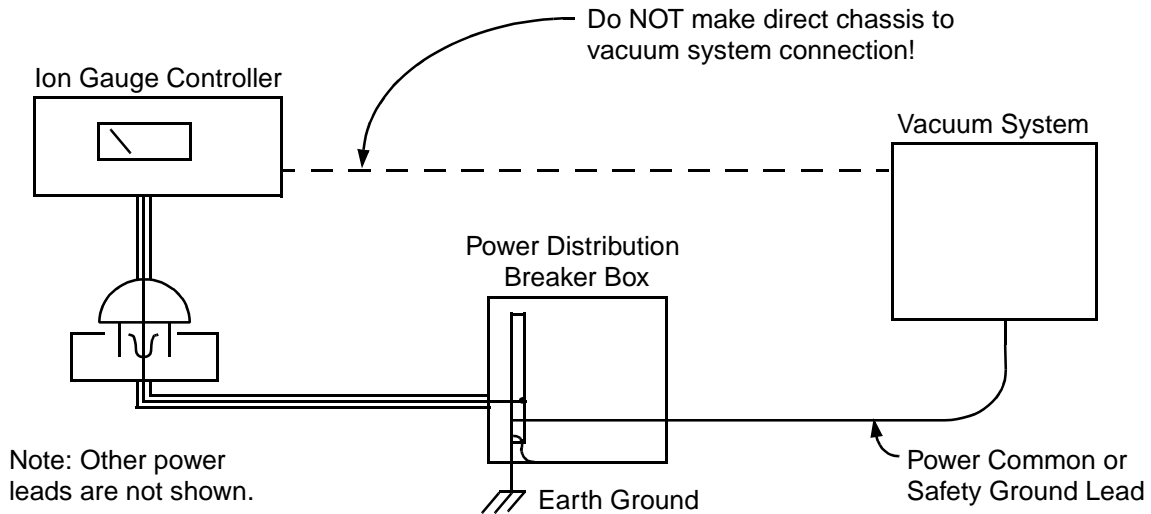


Figure 1 Electrical Connections

An independent agency has determined that all vacuum chambers, regardless of manufacture, can possibly become charged to lethal voltage levels under certain conditions if they are not grounded with a quality, common ground with the controller of their ionization tube.

WARNING



After each maintenance/service procedure and before operating the controller and vacuum system, verify the integrity of the ground of both units.

Equipment utilizing these controls should be designed to prevent personal contact with high voltages.

Always break the primary circuit when direct access to the control unit is required.

Operation

The 563 Ionization Gauge Tube can be used with a variety of commercial controls which are equipped for resistance degassing.



It is particularly important to determine if a control will exceed any of the maximum ratings noted in this manual.

Sensitivity

The sensitivity of an ionization gauge is defined as the ion current per unit of pressure at a specified grid (electron) current. Since all ionization gauges are sensitive to the type and pressure of a gas, the gas composition must be known to correctly establish the pressure reading. Table 1 gives the 563 Ionization Gauge Tube sensitivities for helium, nitrogen and air.

Table 1 Gas Sensitivity Parameters

Gas	Gauge Constant S $S = i^+ / i^- \times 1/P$	Emission Setting for 0.01 Amp/Torr sensitivity (10 uA/micron)
Helium	1.5 per Torr	6.7 mA
Nitrogen	10 per Torr	1.0 mA
Air	10 per Torr	1.0 mA

where:

- S** = Sensitivity, Torr⁻¹
- i⁺** = Ion current, Amps
- i⁻** = Grid (electron) current, Amps
- P** = Partial pressure of the gas, Torr

Specifications

Table 2 Vacuum Specifications

Range	10^{-3} to 2×10^{-10} Torr*
Maximum Operating Pressure	1×10^{-3} Torr
Pumping Speed, ionic	0.06 liters/sec - N ₂ (1 mA)
X-ray Limit*	2×10^{-10} Torr (N ₂) approximately
Initial Calibration	McLeod Standard**

* See Section "X-ray Limit" on page 8

** Calibration of production standards only.

Table 3 Physical Specifications

Envelope	Nonex (all models)
Grid	Tungsten "Non-Sag", 0.025 inch diameter
Filament	Hairpin thoria-coated iridium
Collector	Tungsten, 0.010 inch diameter
Shield Coating	Platinum, internally connected to filament
Base Leads	Soft nickel, 0.060 inch diameter
Collector Lead	Soft nickel, 0.040 inch diameter
Shipping Weight	3 lbs. (6.6 kg)
Internal Volume	220 cc (not including tubulation)

563 Ionization Gauge Tube

Table 4 Operating Ratings

Collector	0 VDC (ground)*
Shield	Internally connected to filament
Grid	+180 VDC
Filament	+30 VDC
Filament Voltage	4.0 VAC
Filament Current	3.5 Amps AC (1 mA grid current)
Filament Voltage, absolute maximum	6.0 VAC
Filament Current, absolute maximum	6.0 Amps
Filament Temperature, absolute maximum	1400 °C

* Collector operated at ground potential through electrometer circuit to reduce noise pickup and leakage currents.

Table 5 Degassing Ratings

Resistance Heating*		Electron Bombardment	
Grid Voltage	6.3 VAC	Grid to Filament Voltage	+700 VDC
Grid Current	8.7 Amps	Grid Current	100 mA DC
Grid Temperature	1200 °C	Grid Temperature	1200 °C max
		Filament Temperature	1400 °C

* All values in the Resistance Heating column are absolute maximum ratings. If these values are exceeded, sagging of the grid can occur.

Mounting

The 563 Ionization Gauge Tube mounts in either a vertical or horizontal position and may be operated and degassed in either position for prolonged periods. An outline drawing of the tube is in Figure 2.

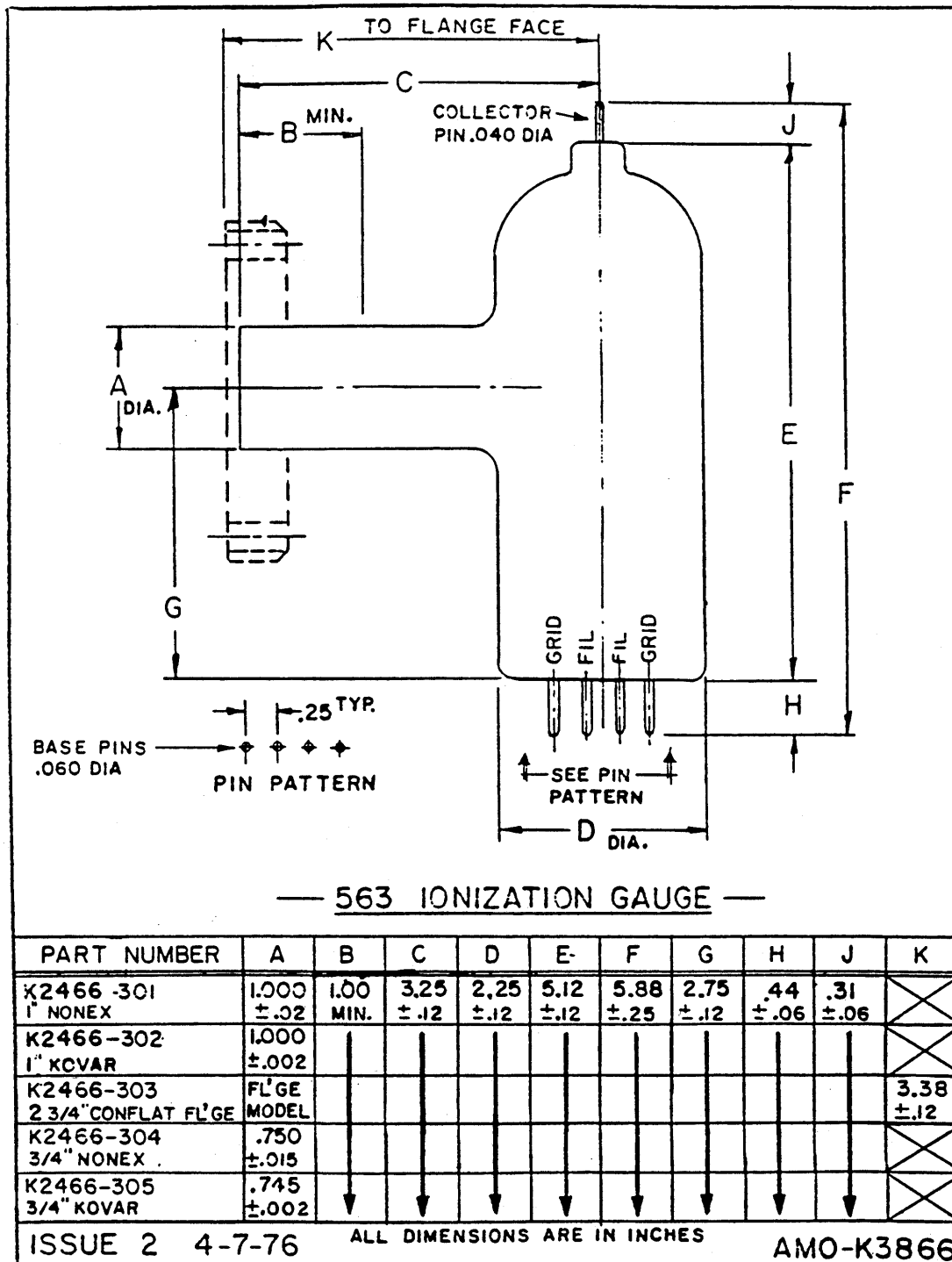


Figure 2 563 Ionization Gauge Tube Outline Drawing

Use of The Gauge Tube

The 563 Ionization Gauge Tube may be used to measure equivalent pressures of gases from 10^{-3} Torr down to the X-ray limit of the gauge. All ionization gauges are composition-sensitive and pressure readings derived from these gauges can only be expressed in terms of equivalent pressure of one gas, usually nitrogen.

Because Bayard-Alpert gauges contain hot filaments, degassing of residual gases in the glass and metal parts of the gauge will occur from the time it is first turned on. It is, therefore, very important to keep the gauge as free from contaminating vapors as possible (diffusion pump fluids, mercury and water vapor, etc.). High temperature baking of the gauge should not be attempted at pressures above 10^{-5} to avoid oxidation of the gauge elements which can make it very difficult to attain ultrahigh vacuum pressure readings. Conversely, ionization gauges cannot be expected to reach low ultimate pressure if the gauge and its connection are not thoroughly degassed.

Degassing

Degassing of the metal parts and glass walls of the gauge can be done by direct resistance heating of the grid or by electron bombardment of the grid by electrons supplied from the filaments of the gauge. In either case, it is important to use the values of grid voltage and current recommended to prevent damage to the gauge (see Table 5 on page 5). The standard Varian controls provide resistance degassing. In general, electron bombardment will result in a faster cleanup of a Bayard-Alpert gauge, but it must be carefully employed and watched to be efficient and non-hazardous to the gauge. Prolonged degassing at pressures above 10^{-7} Torr is usually unnecessary and of little value since the time to re-absorb common gases at this pressure is very short.

Degassing the gauge by resistance heating will typically require one-half hour or more in the 10^{-9} Torr range depending on cleanliness of the vacuum system and its past history. Degassing should proceed until the pressure during degas has reached a peak, then dropped asymptotically to a lower equilibrium pressure. Allowing the gauge to degas overnight will usually insure this if pressures below 10^{-8} Torr are to be measured.

X-ray Limit

In ionization gauges, X-rays generated by the ionizing electrons hitting the grid produce a photoelectric emission at the ion collector. This causes a lower limit of pressure readings known as the X-ray limit (about 2×10^{-10} Torr in the 563 tube). Degassing effects of the gauges are often mistaken for the X-ray limit. Only careful investigation can determine which is the real cause. A second assumption often made is that the X-ray limit is always constant. Among other things, this limit depends on the photoelectric efficiency of the collector surface which, in turn, depends on the amount and type of absorbed gas, etc.

Therefore, as the typical value of X-ray limit is approached (2×10^{-10} Torr nitrogen equivalent), readings should not be considered unequivocal.

Mounting and Flange Availability

Table 6 Ordering Information

Varian Part No.	Description
K2466301	1 inch Nonex
K2466302	1 inch Kovar
K2466303	2 ³¹ / ₄₄ inch Kovar ConFlat® Flange
K2466304	³¹ / ₄₄ inch Nonex
K2661305	³¹ / ₄₄ inch Kovar
K2471311	NW40 KF Flange

This page intentionally left blank.

Gas Correction Factors

Gas correction factor tables are only reproduced for the convenience of the user and do not imply that use with other gases will be safe with hot filament gauge controllers. Table 7 lists relative gauge gas correction factors for various gases. The values are derived by empirical methods substantiated by measurements reported in literature. This table was compiled and published by Robert L. Summers of Lewis Research Center, NASA Technical Note TND-5285, National Aeronautics and Space Administration, Washington, DC, June 1969.

Table 7 Gas Correction Factors

Substance	Formula	Relative Ionization Gauge Gas Correction Factor	Substance	Formula	Relative Ionization Gauge Gas Correction Factor
Acetaldehyde	C ₂ H ₄ O	2.6	Carbon Disulfide	CS ₂	5.0
Acetone	(CH ₃) ₂ CO	3.6			4.7
		4.0			4.8
		3.6	Carbon Monoxide	CO	1.05
Acetylene	C ₂ H ₂	1.9			1.05
		2.0			1.1
Air		1.0	Carbon Tetrachloride	CCl ₄	6.0
		0.98			6.3
Ammonia	NH ₃	1.3	Cesium	Cs	4.3
		1.2			2.0
		1.3			4.8
Amylene:	ISO·C ₅ H ₁₀	5.9	Chlorine	Cl ₂	0.68
		5.8			2.6
Argon	Ar	1.3			1.6
		1.1	Chlorobenzene	C ₆ H ₅ Cl	7.0
		1.2	Chloroethane	C ₂ H ₅ Cl	4.0
		0.9	Chloroform	CHCl ₃	4.7
Benzene	C ₆ H ₆	5.9			4.8
		5.8			4.8
		5.7	Chloromethane	CH ₃ Cl	2.6
		5.9			3.2
		6.0			3.1
Benzoic Acid	C ₆ H ₅ COOH	5.5	Cyanogen	(CN) ₂	2.8
Bromine	Br	3.8			3.6
Bromomethane	CH ₃ Br	3.7			2.7
Butane:	n·	4.9	Cyclohexylene	C ₆ H ₁₂	7.9
		4.7			6.4
	ISO·	4.6	Deuterium	D ₂	0.35
		4.9			0.38
Cadmium	Cd	2.3	Dichlorodifluoromethane	CCl ₂ F ₂	2.7
		3.4			4.1
Carbon Dioxide	CO ₂	1.4	Dichloromethane	CH ₂ Cl ₂	3.7
		1.4	Dinitrobenzene	C ₆ H ₄ (NO ₂) ₂	7.8
		1.5	o·		7.8
		1.5	m·		7.6
		1.4	p·		
			Ethane	C ₂ H ₆	2.6
					2.8
					2.5
			Ethanol	C ₂ H ₅ OH	3.6
					2.9
			Ethyl Acetate	CH ₃ COOC ₂ H ₅	5.0

563 Ionization Gauge Tube

Table 7 Gas Correction Factors, continued

Substance	Formula	Relative Ionization Gauge Gas Correction F5tor	Substance	Formula	Relative Ionization Gauge Gas Correction Factor
Ethyl ether	(C ₂ H ₅) ₂ O	5.1	Naphthalene	C ₁₀ H ₈	9.7
Ethylene	C ₂ H ₄	2.3	Neon	Ne	0.30
		2.4			0.31
		2.2	Nitrobenzene	C ₆ H ₅ NO ₂	7.2
		2.2 to 2.5	Nitrogen	N ₂	1.0
Ethylene oxide	(CH ₂) ₂ O	2.5	Nitrotoluene (o, m, p)	C ₆ H ₄ CH ₃ NO ₂	8.5
Helium	He	0.18	Nitric Oxide	NO	1.3
		0.15			1.2
		0.13			1.0
		0.12	Nitrous Oxide	N ₂ O	1.5
Heptane	C ₇ H ₁₆	8.6			1.7
Hexadiene:					1.7
1.5-	1.5-C ₅ H ₁₀	6.4			1.3 to 2.1
cyclo-	CY-C ₆ H ₁₀	6.0	Oxygen	O ₂	1.0
Hexane	C ₆ H ₁₄	6.6			1.1
Hexene:					0.9
1-	1-C ₆ H ₁₂	5.9			0.9
cyclo	CY-C ₆ H ₁₀	6.4	Pentane		
Hydrogen	H ₂	0.46	n-	n-C ₅ H ₁₂	6.2
		0.38			6.0
		0.41			5.7
		0.45	ISO-	ISO-C ₅ H ₁₂	6.0
		0.44	neo-	(CH ₃) ₄ C	5.7
Hydrogen Bromide	HBr	2.0	Phenol	C ₆ H ₅ OH	6.2
Hydrogen Chloride	HCl	1.5	Phosphine	PH ₃	2.6
		1.6	Potassium	K	3.6
		2.0	Propane	C ₃ H ₈	4.2
		1.5			3.7
Hydrogen Cyanide	HCN	1.5			3.7 to 3.9
		1.6			3.6
Hydrogen Fluoride	HF	1.4	Propene oxide	C ₃ H ₆ O	3.9
Hydrogen Iodide	HI	3.1	Propene:		
Hydrogen Sulfide	H ₂ S	2.2	n-	n-C ₃ H ₆	3.3
		2.2			3.2 to 3.7
		2.3	cyclo-	cy-C ₃ H ₆	3.6
		2.1	Rubidium	Rb	4.3
Iodine	I ₂	5.4	Silver perchlorate	AgClO ₄	3.6
Iodomethane	CH ₃ I	4.2	Sodium	Na	3.0
Isoamyl Alcohol	C ₅ H ₁₁ OH	2.9	Stannic iodide	SnI ₄	6.7
Isobutylene	C ₄ H ₈	3.6	Sulphur Dioxide	SO ₂	2.1
Krypton	Kr	1.9			2.3
		1.7	Sulphur Hexafluoride	SF ₆	2.3
		1.7			2.8
Lithium	Li	1.9	Toluene	C ₆ H ₅ CH ₃	6.8
Mercury	Hg	3.6	Trinitrobenzene	C ₆ H ₃ (NO ₂) ₃	9.0
Methane	CH ₄	1.4	Water	H ₂ O	1.1
		1.5			1.0
		1.6			0.8
		1.4 to 1.8	Xenon	Xe	2.9
		1.5			2.2
		1.5			2.4
Methanol	CH ₃ OH	1.8	Xylene:		
		1.9	o-	o-C ₆ H ₄ (CH ₃) ₂	7.8
Methyl Acetate	CH ₃ COOCH ₃	4.0	p-	p-C ₆ H ₄ (CH ₃) ₂	7.9
Mythyl ether	(CH ₃) ₂ O	3.0			
		3.0			

This page intentionally left blank.



Request for Return Health and Safety Certification



1. Return authorization numbers (RA#) **will not** be issued for any product until this Certificate is completed and returned to a Varian Customer Service Representative.
2. Pack goods appropriately and drain all oil from rotary vane and diffusion pumps (for exchanges please use the packing material from the replacement unit), making sure shipment documentation and package label clearly shows assigned Return Authorization Number (RA#). VVT cannot accept any return without such reference.
3. Return product(s) to the nearest location:

North and South America

Varian Vacuum Technologies
121 Hartwell Ave.
Lexington, MA 02421
Fax: (781) 860-9252

Europe and Middle East

Varian S.p.A.
Via F.lli Varian, 54
10040 Leini (TO) – ITALY
Fax: (39) 011 997 9350

Asia and ROW

Varian Vacuum Technologies
Local Office

For a complete list of phone/fax numbers see www.varianinc.com/vacuum

4. If a product is received at Varian in a contaminated condition, **the customer is held responsible** for all costs incurred to ensure the safe handling of the product, and **is liable** for any harm or injury to Varian employees occurring as a result of exposure to toxic or hazardous materials present in the product.

CUSTOMER INFORMATION

Company name:	
Contact person: Name:	Tel:
Fax:	E-mail:
Ship Method: Shipping Collect #: P.O.#:	
<u>Europe only:</u> VAT Reg. Number:	<u>USA only:</u> <input type="checkbox"/> Taxable <input type="checkbox"/> Non-taxable
Customer Ship To:	Customer Bill To:
.....
.....

PRODUCT IDENTIFICATION

Product Description	Varian P/N	Varian S/N

TYPE OF RETURN (check appropriate box)

<input type="checkbox"/> Paid Exchange	<input type="checkbox"/> Paid Repair	<input type="checkbox"/> Warranty Exchange	<input type="checkbox"/> Warranty Repair	<input type="checkbox"/> Loaner Return
<input type="checkbox"/> Credit	<input type="checkbox"/> Shipping Error	<input type="checkbox"/> Evaluation Return	<input type="checkbox"/> Calibration	<input type="checkbox"/> Other

HEALTH and SAFETY CERTIFICATION

VARIAN VACUUM TECHNOLOGIES CANNOT ACCEPT ANY BIOLOGICAL HAZARDS, RADIOACTIVE MATERIAL, ORGANIC METALS, OR MERCURY AT ITS FACILITY. CHECK ONE OF THE FOLLOWING:

- ☐ I confirm that the above product(s) has (have) **NOT** pumped or been exposed to any toxic or dangerous materials in a quantity harmful for human contact.
- ☐ I declare that the above product(s) has (have) pumped or been exposed to the following toxic or dangerous materials in a quantity harmful for human contact (Must be filled in):

.....

Print Name: Signature: Date:/...../.....

PLEASE FILL IN THE FAILURE REPORT SECTION ON THE NEXT PAGE

Do not write below this line

Notification (RA)#: Customer ID#: Equipment #:

FAILURE REPORT (Please describe in detail the nature of the malfunction to assist us in performing failure analysis):

TURBO PUMPS and TURBOCONTROLLERS

CLAIMED DEFECT <input type="checkbox"/> Does not start <input type="checkbox"/> Noise <input type="checkbox"/> Does not spin freely <input type="checkbox"/> Vibrations <input type="checkbox"/> Does not reach full speed <input type="checkbox"/> Leak <input type="checkbox"/> Mechanical Contact <input type="checkbox"/> Overtemperature <input type="checkbox"/> Cooling defective <input type="checkbox"/> Clogging	POSITION <input type="checkbox"/> Vertical <input type="checkbox"/> Horizontal <input type="checkbox"/> Upside-down <input type="checkbox"/> Other:	PARAMETERS Power: Rotational Speed: Current: Inlet Pressure: Temp 1: Foreline Pressure: Temp 2: Purge flow: <hr/> OPERATION TIME:
Describe Failure :		
TURBOCONTROLLER ERROR MESSAGE:		

ION PUMPS/CONTROLLERS

<input type="checkbox"/> Bad feedthrough <input type="checkbox"/> Poor vacuum <input type="checkbox"/> Vacuum leak <input type="checkbox"/> High voltage problem <input type="checkbox"/> Error code on display <input type="checkbox"/> Other	Describe failure: <hr/> Customer application:
--	--

VALVES/COMPONENTS

<input type="checkbox"/> Main seal leak <input type="checkbox"/> Bellows leak <input type="checkbox"/> Solenoid failure <input type="checkbox"/> Damaged flange <input type="checkbox"/> Damaged sealing area <input type="checkbox"/> Other	Describe failure: <hr/> Customer application:
--	--

LEAK DETECTORS

<input type="checkbox"/> Cannot calibrate <input type="checkbox"/> No zero/high background <input type="checkbox"/> Vacuum system unstable <input type="checkbox"/> Cannot reach test mode <input type="checkbox"/> Failed to start <input type="checkbox"/> Other	Describe failure: <hr/> Customer application:
--	--

INSTRUMENTS

<input type="checkbox"/> Gauge tube not working <input type="checkbox"/> Display problem <input type="checkbox"/> Communication failure <input type="checkbox"/> Degas not working <input type="checkbox"/> Error code on display <input type="checkbox"/> Other	Describe failure: <hr/> Customer application:
--	--

ALL OTHER VARIAN PUMPS

<input type="checkbox"/> Pump doesn't start <input type="checkbox"/> Noisy pump (describe) <input type="checkbox"/> Doesn't reach vacuum <input type="checkbox"/> Overtemperature <input type="checkbox"/> Pump seized <input type="checkbox"/> Other	Describe failure: <hr/> Customer application:
---	--

DIFFUSION PUMPS

<input type="checkbox"/> Heater failure <input type="checkbox"/> Electrical problem <input type="checkbox"/> Doesn't reach vacuum <input type="checkbox"/> Cooling coil damage <input type="checkbox"/> Vacuum leak <input type="checkbox"/> Other	Describe failure: <hr/> Customer application:
--	--

VPD Service Operation

Returned Material Report

This report must accompany all products returned for repair, replacement, or warranty evaluation. Full information regarding reasons for return of the product will expedite repair or adjustment. Please fill in all blanks below and furnish any other information which will help identify the nature and cause of failure.

Reason for Return (check appropriate box)

- | | | | |
|--|---|--|---------------------------------|
| <input type="checkbox"/> Paid Repair | <input type="checkbox"/> Advance Exchange | <input type="checkbox"/> Shipping Error | <input type="checkbox"/> Credit |
| <input type="checkbox"/> Warranty Evaluation | <input type="checkbox"/> Loaner Return | <input type="checkbox"/> Shipping Damage | |

Product Information (use separate forms if more than one model no.)

Varian Model No. _____ Serial No. _____ Quantity _____
Part Description _____

Purchase Information (if product is being returned for warranty evaluation, show your original purchase order number and date purchased)

Varian Sales Order No. (if available) _____ Machine # _____

Original Purchase Order No. _____ Purchase Order Date _____

Company Name _____ Contact _____

Address _____

City _____ State _____ Zip _____

Telephone _____

Failure Report (describe in detail suspected cause or nature of malfunction)

Returned Products

All products returned to Varian/VPD Service Operation for warranty evaluation must be sent **prepaid** and customer must comply with the **warranty replacement and adjustment** provision set forth in the warranty.

Ship directly to: Varian Vacuum Products
Vacuum Products Service Center
121 Hartwell Avenue
Lexington, MA 02421

All products sold by Varian and returned by customer are subject to Varian Vacuum Products standard terms and conditions of sale including, but not limited to, the warranty and damages and liability provisions set forth in the warranty.

Sales and Service Offices

Argentina

Varian Argentina Ltd.

Sucursal Argentina
Av. Ricardo Balbin 2316
1428 Buenos Aires
Argentina
Tel: (54) 1 783 5306
Fax: (54) 1 786 5172

Australia

Varian Australia Pty Ltd.

679-701 Springvale Road
Mulgrave, Victoria ZZ 3170
Australia
Tel: (61) 395607133
Fax: (61) 395607950

Benelux

Varian Vacuum Technologies

Rijksstraatweg 269 H,
3956 CP Leersum
The Netherlands
Tel: (31) 343 469910
Fax: (31) 343 469961

Brazil

Varian Industria e Comercio Ltda.

Avenida Dr. Cardoso de Mello 1644
Vila Olimpia
Sao Paulo 04548 005
Brazil
Tel: (55) 11 3845 0444
Fax: (55) 11 3845 9350

Canada

Central coordination through:

Varian Vacuum Technologies
121 Hartwell Avenue
Lexington, MA 02421
USA
Tel: (781) 861 7200
Fax: (781) 860 5437
Toll Free: (800) 882 7426

China

Varian Technologies - Beijing

Room 1201, Jinyu Mansion
No. 129A, Xuanwumen Xidajie
Xicheng District
Beijing 1000031 P.R. China
Tel: (86) 10 6641 1530
Fax: (86) 10 6641 1534

France and Wallonie

Varian s.a.

7 avenue des Tropiques
Z.A. de Courtaboeuf – B.P. 12
Les Ulis cedex (Orsay) 91941
France
Tel: (33) 1 69 86 38 13
Fax: (33) 1 69 28 23 08

Germany and Austria

Varian Deutschland GmbH

Alsfelder Strasse 6
Postfach 11 14 35
64289 Darmstadt
Germany
Tel: (49) 6151 703 353
Fax: (49) 6151 703 302

India

Varian India PVT LTD

101-108, 1st Floor
1010 Competent House
7, Nangal Raya Business Centre
New Delhi 110 046
India
Tel: (91) 11 5548444
Fax: (91) 11 5548445

Italy

Varian Vacuum Technologies

Via F.lli Varian, 54
10040 Leini, (Torino)
Italy
Tel: (39) 011 997 9111
Fax: (39) 011 997 9350

Japan

Varian Vacuum Technologies

Sumitomo Shibaura Building, 8th Floor
4-16-36 Shibaura
Minato-ku, Tokyo 108
Japan
Tel: (81) 3 5232 1253
Fax: (81) 3 5232 1263

Korea

Varian Technologies Korea, Ltd.

Shinsa 2nd Bldg. 2F
966-5 Daechi-dong
Kangnam-gu, Seoul
Korea 135-280
Tel: (82) 2 3452 2452
Fax: (82) 2 3452 2451

Mexico

Varian S.A.

Concepcion Beistegui No 109
Col Del Valle
C.P. 03100
Mexico, D.F.
Tel: (52) 5 523 9465
Fax: (52) 5 523 9472

Taiwan

Varian Technologies Asia Ltd.

18F-13 No.79, Hsin Tai Wu Road
Sec. 1, Hsi Chih
Taipei Hsien
Taiwan, R.O.C.
Tel: (886) 2 2698 9555
Fax: (886) 2 2698 9678

UK and Ireland

Varian Ltd.

28 Manor Road
Walton-On-Thames
Surrey KT 12 2QF
England
Tel: (44) 1932 89 8000
Fax: (44) 1932 22 8769

United States

Varian Vacuum Technologies

121 Hartwell Avenue
Lexington, MA 02421
USA
Tel: (781) 861 7200
Fax: (781) 860 5437
Toll Free: (800) 882 7426

Other Countries

Varian Vacuum Technologies

Via F.lli Varian, 54
10040 Leini, (Torino)
Italy
Tel: (39) 011 997 9111
Fax: (39) 011 997 9350

Internet Users:

Customer Service & Technical Support:

vppl.customer.support@varianinc.com

Worldwide Web Site:

www.varianinc.com/vacuum

Order On-line:

www.evarian.com

Representatives in most countries

